

Model 3970



## Traverse rail with connector

Connect the traverse system to a second room




- Accessing an adjoining room is possible with a traverse connection.
- Allows the motor to move smoothly from one room to another by connecting the movable rail to a second room with a fixed rail.
- A manual version is easy to use with little effort from the caregiver.
- An electric version is available so the entire system can be used independently.
- An integrated safety allows the motor to continue only when both rails are perfectly connected.


Lift & Care systems for people with reduced mobility

[www.handimove.com](http://www.handimove.com)

# Traverse rail with connector

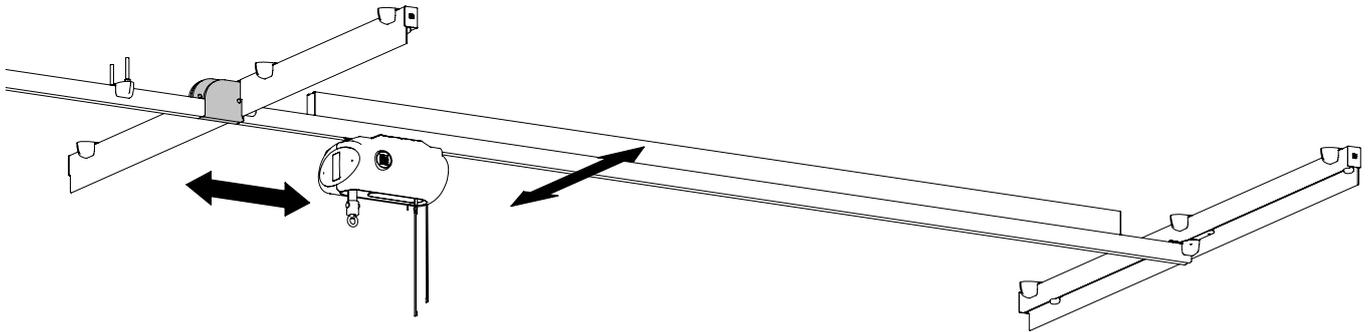


## ● Dimensions (mm)

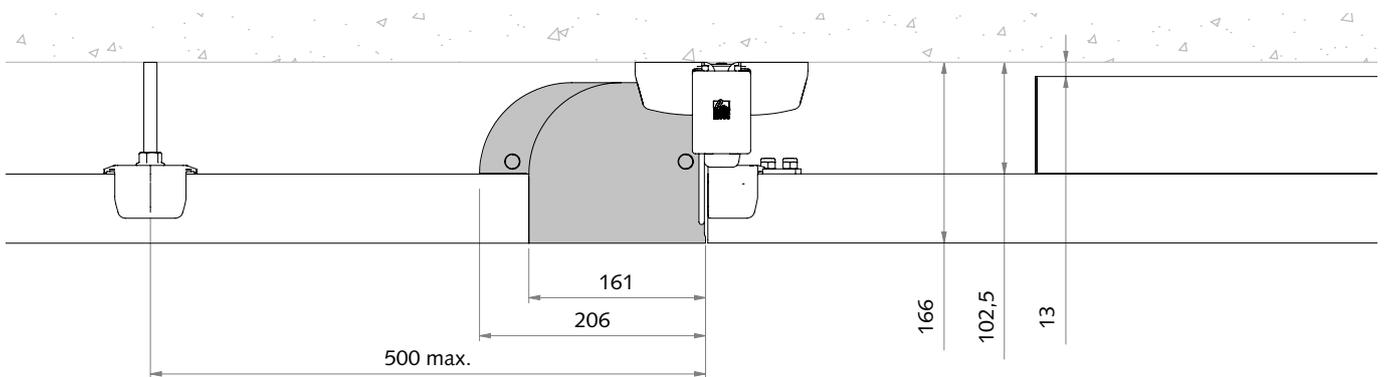

\* The drawing is an example: for a traverse rail smaller or equal to 5 meters.

## ● Technical specifications

Safe working load.....240 kg  
 Material.....Aluminum painted, stainless steel, ELVZ steel, POM, FKM

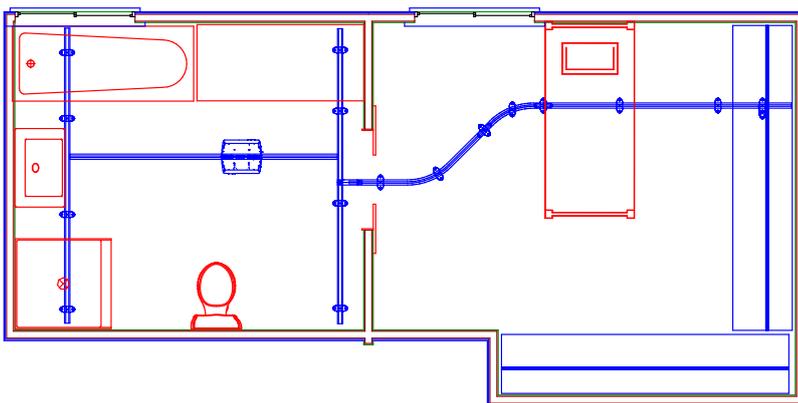

For free demonstration or further information

