

Model 2800-2805

# Ceiling motor

Efficient and fast



- The fastest track motor
- Suitable for intensive use
- Featherlight sideways movement
- To be used with slings and/or the Handi-Move® Body Support
- To be used with classic spreader bar or electric tilting spreader bar
- Motor and rail systems make it possible to access all areas
- Wireless control optional
- Back-up control on motor
- Emergency down and stop within reach

Lift & Care systems for people with reduced mobility


[www.handimove.com](http://www.handimove.com)

# Ceiling motor

## ● Technical specifications

Safe working load .....	200 kg/300 kg (3000 mm motorband: 160 kg)
Batteries .....	2 x 12 V 9.5 Ah/ rechargeable
Standard motor band length .....	2000 mm (option: 3000 mm)
Rail System .....	aluminium, lacquered in white
Smooth stop and go .....	electronic
Cover .....	ABS
Emergency down .....	mechanical
Handcontrol 2800 .....	watertight IP67 electric
Handcontrol 2805 .....	infra-red
Number of lifts with fully charged batteries .....	160
Load and battery indications .....	visual and auditive
Charger .....	input 0.9 A 100-240 Vac, 50-60 Hz Class II
.....	output 24 Vdc Max 1.3 A (9940) Max 1.2 A (2541)

## ● Examples of routes

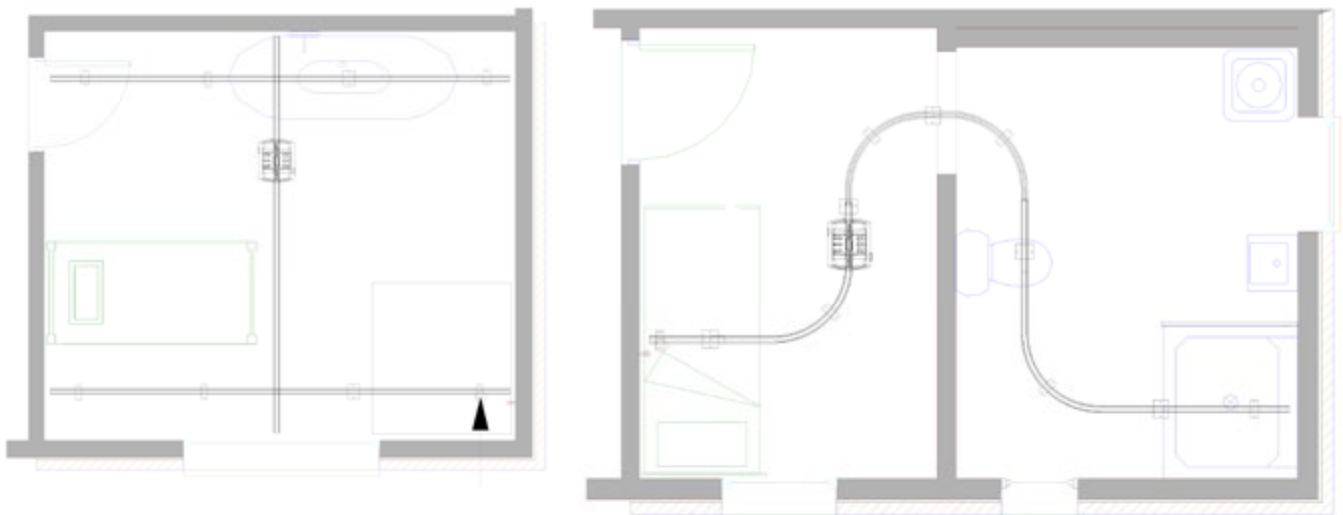

## ● Versions\*

### 2800

Lifting motor (up/down)  
With electric handcontrol

### 2805 IR

Lifting motor (up/down)  
With infra-red remote control

\* also available with lateral motor (2810 L - 2815 L IR)

For free demonstration or further information

