

Model 2840

Portable ceiling hoist

A portable motor for multiple transfers in several rooms:
you can use it everywhere





- Can be used in different rooms with different rail sections (bedroom, bathroom, toilet, shower...).
- The door opening does not have to be modified.
- Light and compact, yet powerful: lifting up to 182 kg.
- Aesthetic: streamlined timeless design.
- Budget-friendly: one motor for several places.
- In combination with the Handi-Move Body Support® or slings.

Lift & Care systems for people with reduced mobility


www.handimove.com

Portable ceiling hoist

● Dimensions

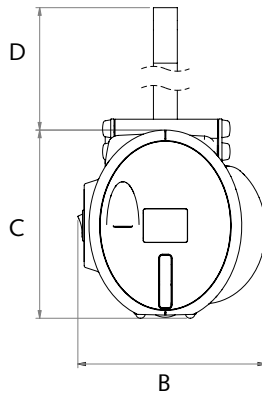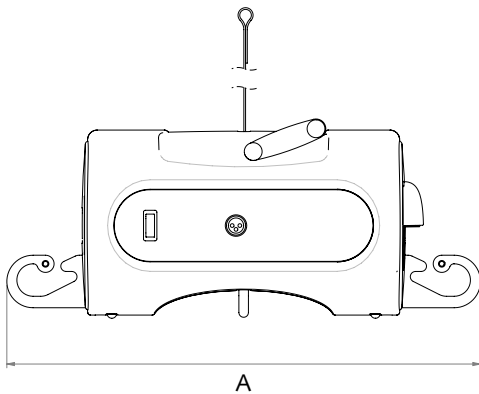

A	480 mm
B	200 mm
C	210 mm
D	1550 mm

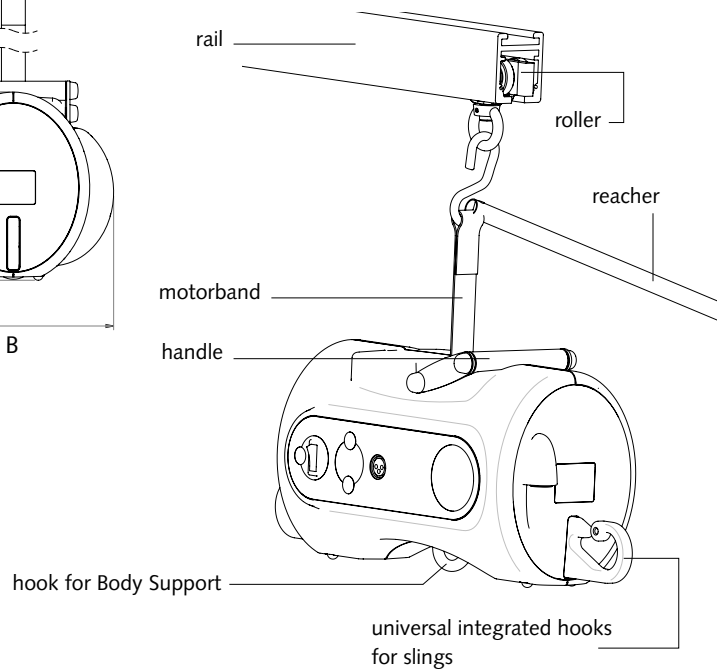

● Technical specifications

Weight	8 kg
Max. Lifting Capacity	182 kg
Batteries	2 x 12 VDC/3.2 Ah.
Charger:	input 220 V AC / output 24 V DC
Number of lifts with full batteries.....	50 with 100 kg / 50 cm
Motor cover	ABS
Emergency down	Electric
Overhead rail.....	Aluminium alloy, white laquered

● With slings or Body Support

Thanks to the integrated hooks on both sides of the motor and the swivel hook placed centrally at the bottom, the motor can be used with all slings or with the Handi-Move Body Support.



Option: Transport trolley

For free demonstration or further information