
INTERNATIONAL

Model 3200

Pool lift

Direct from wheelchair to pool

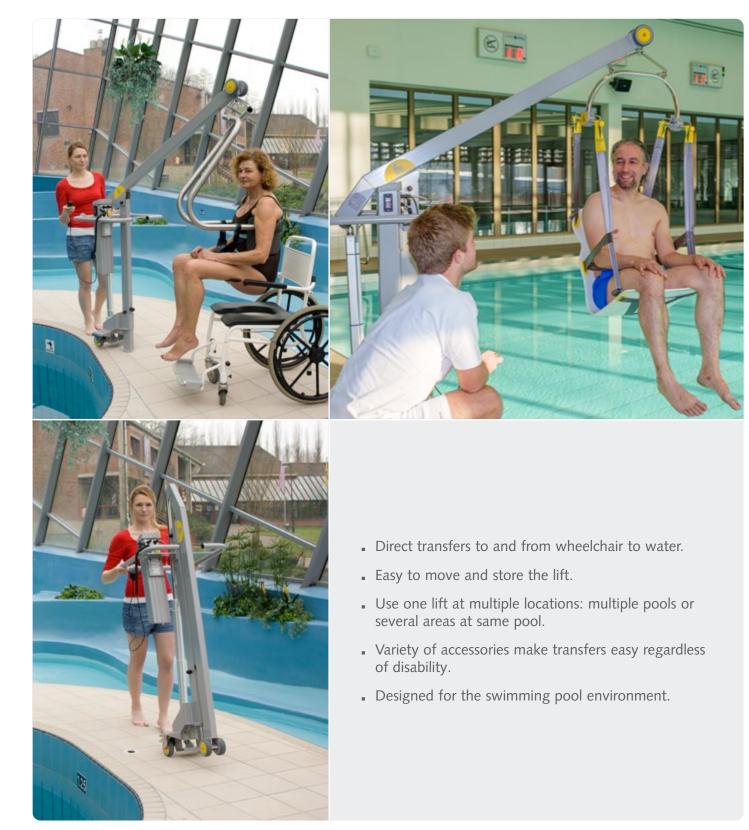

Lift & Care systems for people with reduced mobility

Model 3200



Pool lift

DIRECT TRANSFERS TO AND FROM WHEELCHAIR TO WATER.

The Pool lift brings the user to the water in one smooth motion. This efficient lift moves the user to the pool in one transfer.



The Pool lift has a pivoting range of 360°.

EASY TO MOVE AND STORE THE LIFT.

Operation is simple, with the push of a button, the spindle lowers into the floor socket at the edge of the pool.

When finished, lift the spindle out of the socket, and wheel the lift to storage.

Once removed, seal the socket opening with the plug.

USE ONE LIFT AT MULTIPLE LOCATIONS: MULTIPLE POOLS OR SEVERAL AREAS AT SAME POOL.

The installation of additional floor sockets, and the fact that it is easily movable, allows the lift to be used in different locations


VARIETY OF ACCESSORIES MAKE TRANSFERS EASY REGARDLESS OF DISABILITY.

Customize the transfer with the right accessory.

Based on the need choose a pool lift chair, Handi-Move Body Support, bathing sling or standard sling.

DESIGNED FOR THE SWIMMING POOL ENVIRONMENT.

Made of stainless steel and powered by batteries, the mobile Pool lift is designed for safe, long-term use in the pool environment.

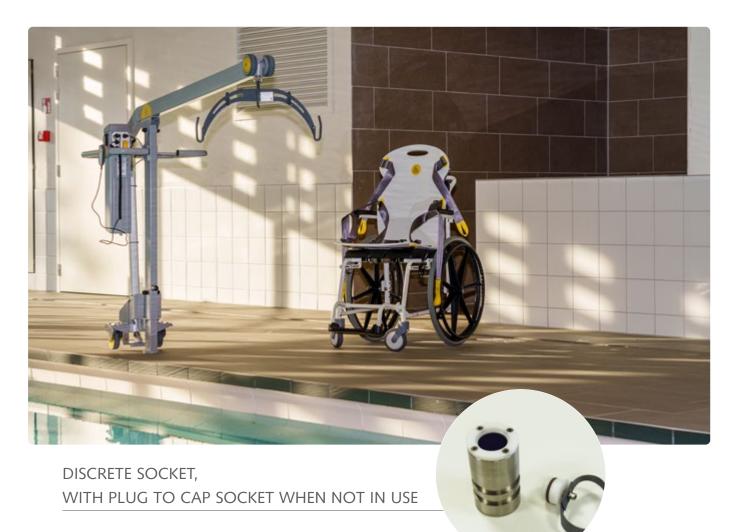

© 31-05-2024 • Handi-Move • v1.4

Model 3200


Pool lift





DIMENSIONS AND SPECIFICATIONS


Dimensions (mm)		
А	1840	
В	1305	
С	540	
D	350	
E	1700	
F	550	
G	500	
Н	310	
I	370	

Specifications	
Safe working load	135 kg
Batteries	2 x 12 Vdc; 8 Ah
Charger	separate
Diameter of casters	100 mm
Construction	stainless steel 316
Hand control	electrical (IP67)

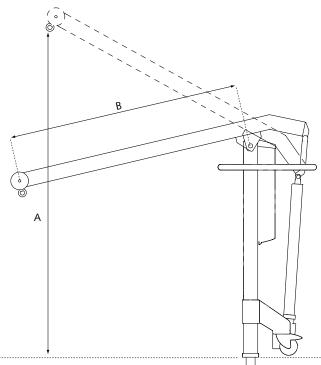

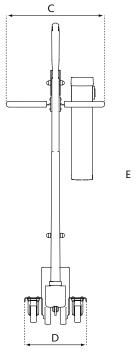

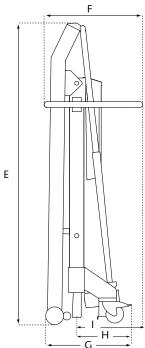

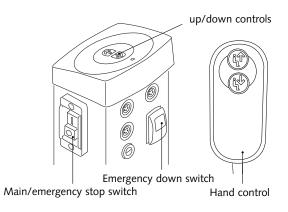


© 31-05-2024 • Handi-Move • v1.4


Model 3200


Pool lift



BODY SUPPORT SYSTEMS

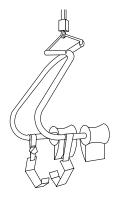

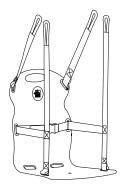

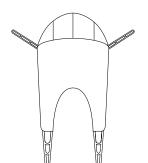

Body Support

XS - S - M • art. 1035 Electrolytically coated (anticorrosive) + Slipproof Thigh Supports • art. 1074



Pool Lift Chair (PVC)

Uni-size • art. 3620 + Spreaderbar • art. 2110 / 2125




Standard Sling with moon shaped head support and padded thigh flaps

XS - S - M - L - XL • art. 3431 + Spreaderbar • art. 2110 / 2125



Bathing Sling (PVC) XS - S - M - L - XL • art. 3464

+ Spreaderbar • art. 2110 / 2125

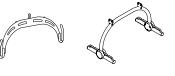

For free demonstration or further information

© 31-05-2024 • Handi-Move • v1.4